

KASON KS-VM-3 Vacuum Mounting Machine



(Picture just for reference)

19 YEARS

Development Base of Advanced Testing Machine Of China

Application

KS-VM-3 Vacuum Mounting Machine is a special equipment for cold mounting of samples. The sample and mounting resin can be degassed by evacuating the sample and mounting resin. In this way, the voids in the sample will be filled with the mounting resin, and there will be no gaps between the sample and the mounting resin. The device is compact, fast and easy to use. 10 samples with a diameter of 30mm can be mounted at the same time.

Technical Parameter

Name	Specification	Remarks
Die specification	φ30mm, soft type	Optional φ40, φ50mm
Vacuum pump	Maximum vacuum: -0.1Mpa, Pumping rate: 7.2 m ³ /h	
Voltage and power	Single phase 220V, 50/60Hz 370w	
Dimension	Vacuum container: 300L*300W * 300H mm Vacuum Pump: 340L* 130W * 270H mm	
Maximum number of mountable samples	>10 PCS	φ30mm
Weight	Vacuum container: 1.15kg Vacuum pump: 7.5kg	

Configuration

Name	Quantity
Host	1 Set
Vacuum pump	1 Set
Vacuum tube	1
High transparent epoxy king sample	1 Set
Mold φ30	1
Manual	1
Disposable plastic cups	10 Pcs

19 YEARS

Professional focused on testing equipment

KASON is established in 2003, owns more than 8000 square meters factory. has a professional sales teams, modern enterprise technology center, scientific and technological research and development team.


Machines passed the European CE authentication, American FDA certificate and ISO 9001.

Products sold to USA, Canada, Australia, Europe, Africa etc, more than 130 countries and supply OEM service for many customers

PROFESSIONAL TEAM

KASON has a professional sales teams, modern enterprise technology center, scientific and technological research and development team.



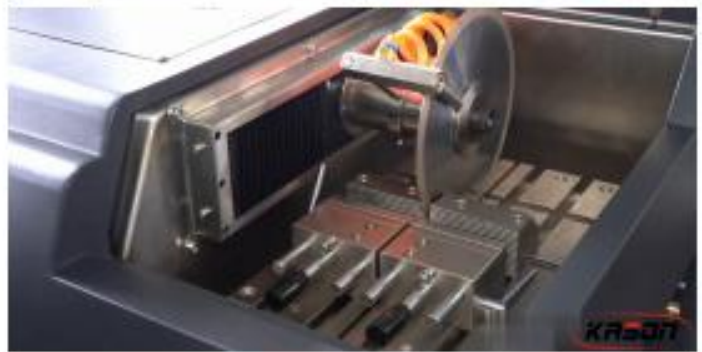

KASON reserves the right to modify the technical and stetics characteristics included in this document, without previous notice